

ORDER OF DRAW IMPORTANT!!!

CONTAMINATION CAN PRODUCE INCORRECT RESULTS


To avoid **EDTA** contamination, samples **MUST** be taken in the following order:

1. Blood culture (if necessary)

2. Serum



3. Serum Gel

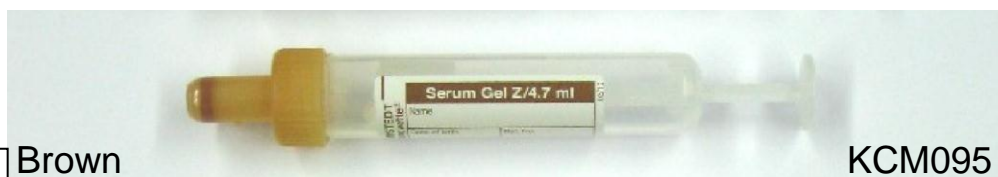

4. Citrate

Fill level is critical to analysis
Must be filled to mark
(↓)



5. Lithium Heparin




6. Potassium EDTA



7. Blood Transfusion EDTA



8. Fluoride EDTA

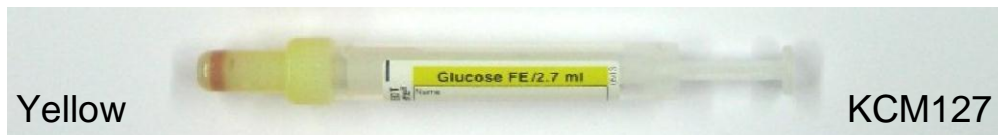

FILL TUBES UP TO THE MARK AND MIX BY INVERSION
DO NOT POUR BLOOD FROM ONE TUBE TO ANOTHER

Document Code:	Owner: J Sullivan
Version Number: 1.0	Author: J Sullivan
Date of Issue:	Authority for Issue: L Bailey Signature: L Bailey